

°CRYO SCIENCE
THE °CRYO ARCTIC™



°CRYO ARCTIC™

The °CRYO Arctic™ is designed to deliver the easiest, safest, and most effective whole body cryotherapy treatments in the industry. Built using state of the art technology, and created to ensure that the client is never in direct contact with nitrogen vapors, the °CRYO Arctic™ matches its futuristic design with the highest safety standards.

The user-friendly and streamlined interface is equipped with a library of °CRYO protocols to ensure a simple yet effective user journey.



°CRYO ARCTIC™ - Features



°CRYO ARCTIC™ - Features Comparison

	Features	°CRYO Arctic™
Treatment Types	Whole Body Cryotherapy (WBC)	✓
	Partial Body Cryotherapy (PBC)	✓
Treatments	Air Flow Control	✓
	Temperature Control	✓
	°CRYO Pro Mode	✓
Safety Features	Air Flow Control	✓
	Temperature Control	✓
	°CRYO Pro Mode	✓
Features	Digital touchscreen to control the °CRYO Arctic	15.6"
	Client facing screen to display time remaining and marketing text	✓
	Window	Automatic
	Troubleshooting with Wi-Fi capability	✓
	High quality speakers to play music via Bluetooth	✓
	Audio communication system between operator and client	✓
	Interior Finish	Dual Tone
Nitrogen Vessel	Average Liquid Nitrogen Consumption	12L per session

°CRYO ARCTIC™ - Technical Specifications




°CRYO Arctic™	
Height	2300 mm / 91 in
Width	1800 mm / 71 in
Depth with door	1000 mm / 41 in
Number of shipping crates	3
Total Weight without crates	410 kg / 904 lbs
Total Weight with crates	750 kg / 1655 lbs
Cooling Medium	Cryogenic Liquid Nitrogen (LIN)
Temperature within chamber	-110 to -140°C / -166°F to -266°F
Electrical requirements	230 V, 50Hz/110V, 60Hz


°CRYO ARCTIC™ WHOLE BODY CHAMBER COOLING PROCESS

HIGHEST SAFETY STANDARDS - NO DIRECT CONTACT WITH NITROGEN



Unlike other cryotherapy chambers, the °CRYO Arctic™ only exposes the client to breathable cold air, and not to harmful Nitrogen gas. Nitrogen is used only as a cooling medium and does not come into contact with a person inside.

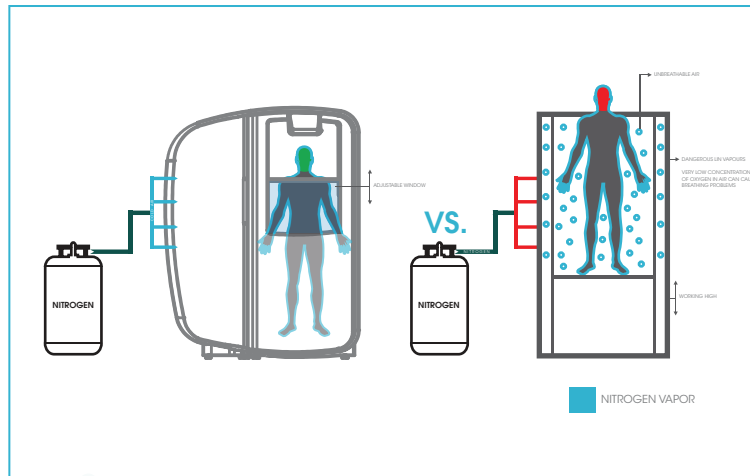


The °CRYO Arctic™ has more space for clients to move around inside the chamber and not feel claustrophobic.

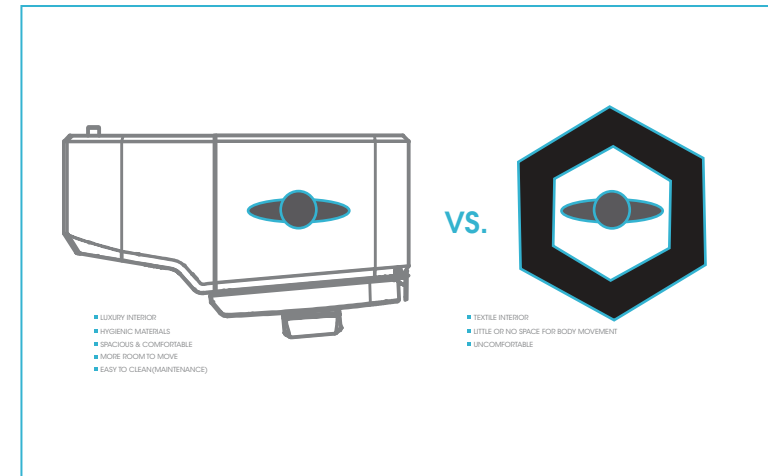

The °CRYO Arctic™ has an adjustable window in case the client would like to have his head above the air vapor.

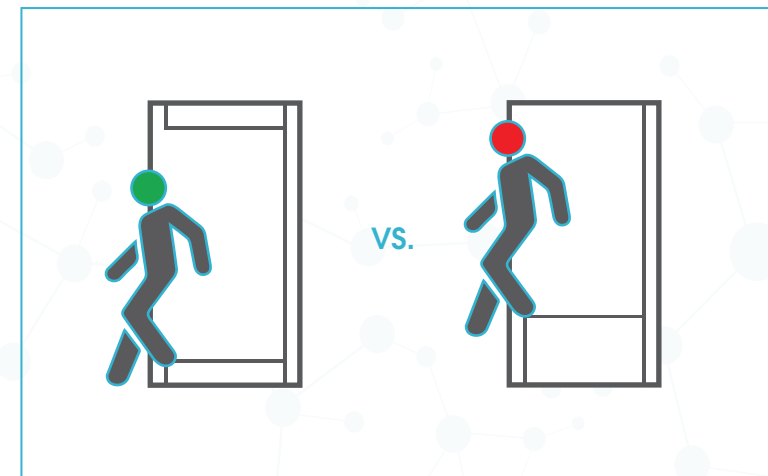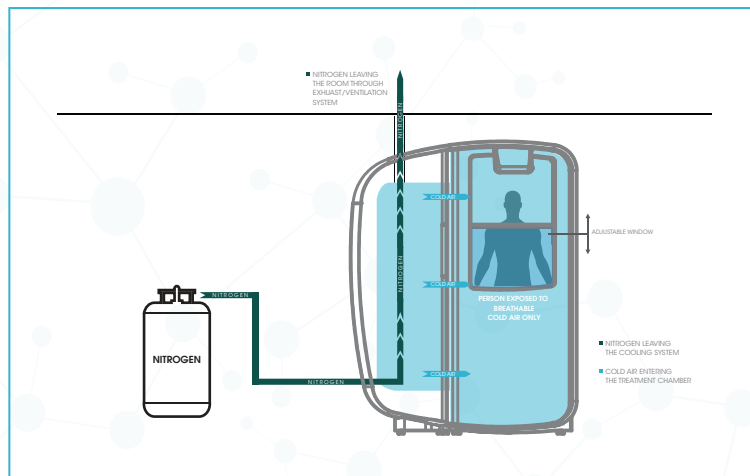

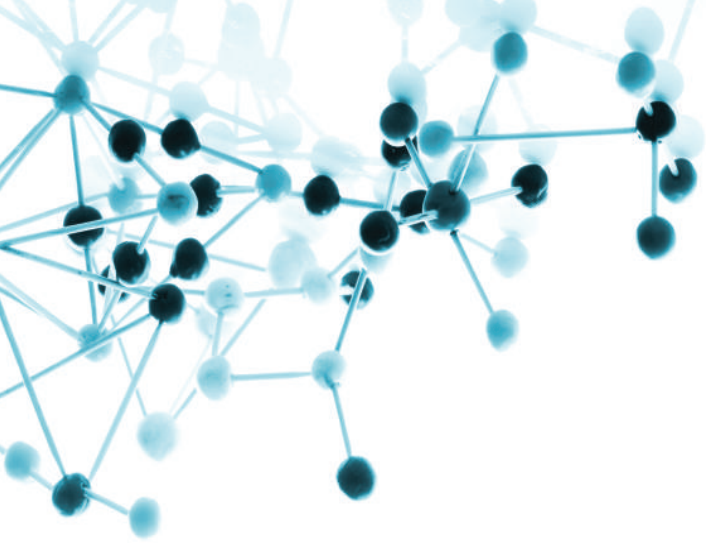

www.cryosc.com

Get in touch with us

contact@cryosc.com

www.cryosc.com

